

## TAURUS ...



#### **Features**

- Allows just one functioning handpiece at a time
- While operating with a handpiece, no accidental movement of the chair
- Emergency halt with a press of any key
- Fiber Optic upgrade available at no extra cost.
- Free from pressing the disk pedal for continuous work with a continuously functioning scaler. Pre-setting option available.

## Control Switch Panel

 Easy control of functions thanks to the direct buttons and intuitive color scheme

LCD Display (Option





### **Doctor Table**

• With an open-type holder allowing left and right rotation, it easily secures convenient positions for instrument usage.





### Soft Chair Movements

• Wide range of chair height - between 16±0.8" and 27±0.8" - and soft chair movements offer the best positions for both patients and operators



## Twin-articulating headrest

 Various angles are enabled on the headrest, providing optimum positions for both adults and children (optional)





## Rotating Armrest

 The armrest ensures zero obstruction to patient entry and exit; also provides comfort tor the patient's arm



## Easy-to-Clean Spittoon

 Designed for a hygienic environment, the spittoon is easily detachable and easy to clean.




# Rotating chair-mounted cuspidor unit

• The compact swiveling body of the device allows the assistant to have better access to the patient

## Specification

- Rated voltage and frequency: AC 110-120V~, 50/60HZ
- Power consumption: 1,000VA
- Seat movement powered by hydraulic motor
- Chair elevation range: 16±0.8" and 27±0.8"
- Backrest inclination: 0°±3° ~ 70°±3°
- Seat length: 2.8"
- Weight(net): 366±11lbs



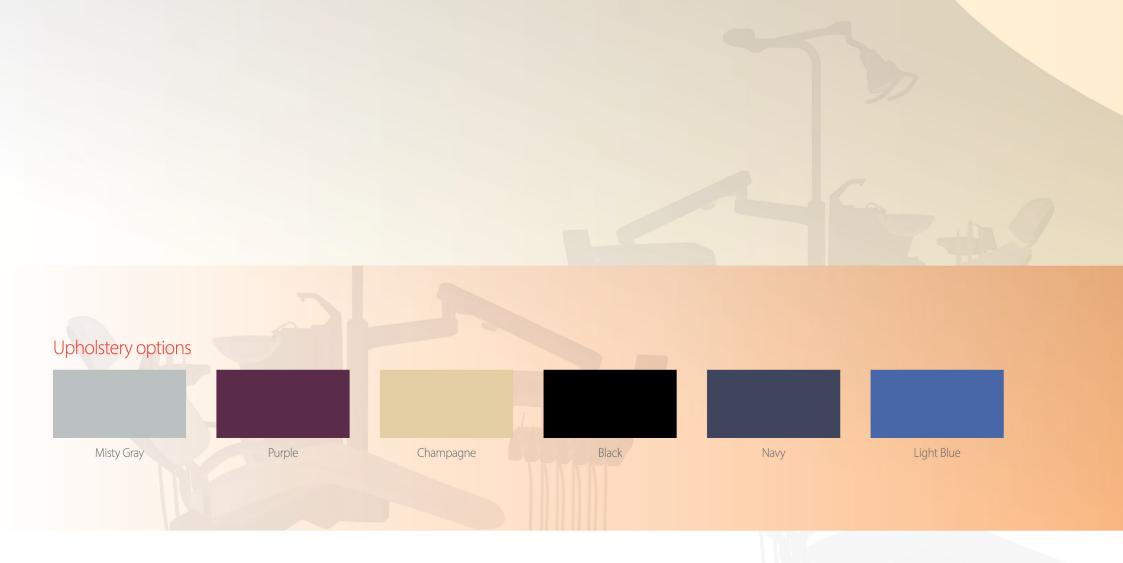

### Manufacturer

#### Shinhung Co., Ltd.

450 CHEONGPA-RO, JUNG-GU, SEOUL, KOREA Tel: 82-2-6366-2041 / Fax: 82-2-6366-2131 www.shinhung.com

E-mail: export@shinhung.co.kr

We reserve the rights to modify the design or specification without prior notice to allow for the improvement.

## **Distributor**

#### SIS-DENTAL INC

6131 ORANGETHORPE AVE SUITE 410, BUENA PARK, CA 90620, USA

TEL: 714-752-6001 www.shinhungimplant.com

E-mail: info@sis-dental.com

CC/1511/2000 Rev.0